

# YOUR OWN PRESCRIPTION SOLUTION

CLEAR VISION ALLOWS YOU TO SEE EVERY TINY DETAIL, SO THAT YOU CAN PERFORM BETTER AND ENJOY THE SPORTS YOU PLAY.
SEEING CORRECTLY REDUCES EYE FATIGUE, INCREASES CONCENTRATION, AND LOWERS YOUR RISK FOR INJURIES. JUST A SMALL
CHANGE IN PRESCRIPTION DRASTICALLY INCREASES YOUR RISK OF LOSING FOCUS. HIGH QUALITY SPORTS EYEWEAR PROTECTS THE
EYES AND PROVIDES FIRST CLASS VISION. FOR ATHLETES WHO REQUIRE PRESCRIPTION LENSES, ADIDAS SPORT EYEWEAR OFFERS
A RANGE OF PRESCRIPTION SOLUTIONS COMPATIBLE WITH THEIR INDIVIDUAL NEEDS AND TYPES OF SPORT.

YOU CAN CHOOSE BETWEEN DIRECT GLAZING, ADAPTER, AND CLIP-IN SOLUTIONS TO ENSURE OPTIMAL PERFORMANCE THROUGH OPTIMAL VISION – YOUR INDIVIDUAL PRESCRIPTION SOLUTION, BUILT FOR PERFORMANCE!



#### **CLIP-IN**

THIS VERSATILE PRESCRIPTION SOLUTION ALLOWS YOU TO QUICKLY AND EASILY CHANGE LENSES.



#### **ADAPTER**

PROVIDES YOU WITH A WIDE RANGE OF PRESCRIPTION OPTIONS.



#### **DIRECT GLAZING**

UNCOMPROMISED VISION IN SPORT. PROVIDES YOU WITH THE LARGEST FIELD OF VISION FOR PERFECT SIGHT FROM EVERY ANGLE.



#### **PERFECT SIGHT** FOR ATHLETES

The prescription concept by adidas Sport eyewear guarantees sportsmen and women clear and undistorted vision. Built to face anything, with all the performance benefits an athlete needs, with prescription included. This means best performance even in extreme situations, which can be adjusted to individual preferences and diverse sports.

#### THE PRESCRIPTION CONCEPT

adidas Sport eyewear guarantees athletes clear, undistorted vision for maximum performance in every situation, even in extreme locations.

#### AS INDIVIDUAL AS YOU

The prescription concept is adjusted to individual needs and different sport types. For athletes who want to achieve maximum performance with optimum vision.

#### SAFE AND ADVANTAGEOUS

Clear vision means less risk of injury thanks to optimum reaction times and better risk assessment. Reduced eye fatigue increases concentration and visual performance, which results in a higher fun factor in sport.

ONLY THOSE WHO CAN FOCUS ON THEIR TARGET WILL ACHIEVE IT.

The tilted position of the lenses in sharply curved eyewear causes optical side-effects.

In order to obtain the right eyewear, the lenses should be corrected. You can download the correction program Silhouette wrap-around calculator from our website (http://b2b.silhouette.com) free of charge. The program integrates all customer data, frame data and lens data into the calculation, in order to obtain the best possible customer-specific correction of all optical aberrations. Please contact your local distributor. The vertical lens size means that bifocal and progressive lenses can be fitted. However, the lens type must be tailored to the wearer's individual requirements.

#### PLEASE NOTE:

- Correct the distance refraction using the Silhouette wrap-around calculator
- Do not correct the addition (= difference between the distance prescription and the near prescription)

If using bifocal or progressive lenses, the large face angle may make the lenses difficult to wear, so such glazing is carried out under your own responsibility.

#### **RX SPARE PARTS**

FK 85-3 box

SHS17 microfiber bag
P 0029 mounting box

RX clip-in rimless

BLS 63 plastic sleeve
Silicon 1,0 shock absorber

RX adapter nylor

NY 2 nylon thread
AP11 metal cover panel

#### RX DIRECT INFRAME GLAZING





## DIRECT INFRAME GLAZING (up to 10-base frames)

A more permanent solution offering perfect vision in all directions, especially with freeform technology, which is exactly calculated at every angle to the individual needs. Various adidas models offer the possibility of direct inframe clazing.






#### **RX ADAPTER NYLOR**

For adapter glazing.

- wide field of vision
- glazing can be done by an optician (like nylor glazing)
- more glazing possibilities with adapter than with direct glazing
- we recommend polycarbonate or HIX lenses



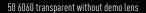



50 6050 black without demo lens


#### RX LENSES

Base 8 compensated due to the high wrap angle (>10^°).



#### **ADAPTATION**

Glazing type: nylor. Glazing with groove depth 0,75 mm, groove width min. 0,75 mm (see diagram). A lens former is available.

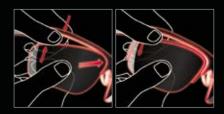

#### MOUNTING

Remove adapter from the frame. Insert RX lens first nasal than temporal in the adapter. Pull the nylor in the groove from frontside.





#### RX ADAPTER FULLRIM

RX optical adapter incl. clear demo lenses for alignment available. Please contact your local distributor.



50 6050 grey without lens extended: 51 6050 grey without lens



50 6060 grey silver mirror without lens extended: 51 6060 grey silver mirror without lens

#### **RX LENSES**

Base 8 compensated due to the high wrap angle (>10^°).



#### **ADAPTATION**

Glazing type: nylor. Glazing with groove depth 0,6 mm, groove width min. 0,6 mm (see diagram). A lens former is available.



#### MOUNTING


adapter + RX lens RX lens inserted into adapter from behind - beginning nasally






#### **RX CLIP-IN RIMLESS**


- glazing parameters approx. -6/+6 HIX lens
- practically invisible
- extensive field of vision
- RX clip-in rimless is protected against external influences by the sun lenses
- most cost effective method for changing sunglass lenses according to different light conditions
- can be used for sport eyewear up to 10-base
- great variety of shapes
- suitable for bifocal and progressive lenses due to the deeper lens shape (a714/a715)

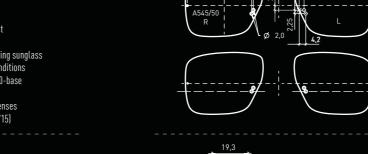




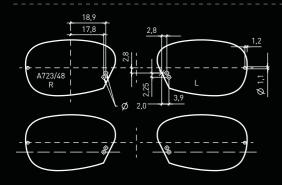




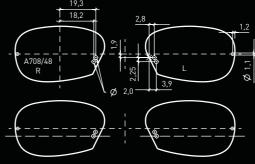



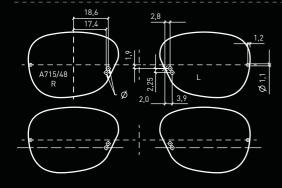



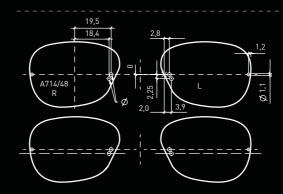



RECOMMENDATION

high index lens 1,6 | dioptries:+/-6 dpt possible | antireflective coating

















a727	
*	39 x 32 mm (921 mm²)
	39 mm
X	26 mm
	10°
Attention: short	distance to eyelashes!



<b>.</b>	
×	46 x 31 mm (1099 mm²)
	47 mm
X	19 mm
	17°



#### RX CLIP-IN RIM

- glazing parameters approx. -10/+8RX clip-in rim is protected against external influences by the sun lenses
- most cost effective method for changing sunglass lenses according to different light conditions
- can be used for sport eyewear up to 10-base
- easiest glazing
- easiest mounting (warming not necessary)
   the most cost-saving solution regarding the lenses
- independent of lens material
- quick and easy lens change



- CR 39/ high index lens 1,6 this depends on dioptries
- dioptries:+/-10 dpt possible
- antireflective coating

a546











	39 x 32 mm (921 mm²)
	39 mm
X	26 mm
	10°

# PRESCRIPTION-OVERVIEW

									+/-
			RX	CLIP-IN RIMLES		RX CLIP		RX ADAPTER	RX DIRECT INFRAME GLAZING
MODEL	NAME	SIZE	CHASSIS	fixed with OPTION 1	demo lens OPTION 2	fixed with demo lens OPTION 1 OPTION 2			IN RAIVIL GLAZING
ad01	ZONYK PRO	L	a515	a545					
ad02	ZONYK PRO	S	a515	a545					
ad03	ZONYK	L	a515	a545					
ad04	ZONYK	S	a515	a545					
ad05	ZONYK AERO PRO	L/S	a515	a545					
ad06	ZONYK AERO	L/S	a515	a545					
ad20	JAYSOR		a515	a714		a779			
ad21	BABOA		a515	a714		a779			
ad22	HORIZOR		a515	a714		a779			
ad80	BACKLAND					a546			
ad81	PROGRESSOR C					a546			
ad82	PROGRESSOR S					a546			
ad83	PROGRESSOR PRO					a546			
ad84	BACKLAND DIRT					a546			
a162	ID2					a747			
a167	EVIL EYE HALFRIM PRO	L	a515 / a516	a714	a715	a731	a779	a530	
a168	EVIL EYE HALFRIM PRO	S	a515 / a516	a714	a715	a731	a779	a531	
a180	EVIL EYE HALFRIM PRO	XS	a515 / a516	a708	a715	a731	a779	a539	
a181	EVIL EYE HALFRIM PRO	L	a515 / a516	a714	a715	a731	a779	a530	
a184	ID2 PRO					a747			
a189	TYCANE PRO	L	a515 / a516	a714	a715	a731	a779	a540	
a190	TYCANE PRO	S .	a515 / a516	a714	a715	a731	a779	a541	

# PRESCRIPTION-OVERVIEW

a191 TYC/ a192 TYC/ a192 TYC/ a193 EVIL	ANE L EYE EVO PRO L EYE EVO PRO SANE PRO OUTDOOR	L a5'	HASSIS C	a714 a714	s DPTION 2 	RX CLIP fixed with d OPTION 1  a731		RX ADAPTER  a540	RX DIRECT INFRAME GLAZING
a191 TYC/ a192 TYC/ a192 TYC/ a193 EVIL	ANE L EYE EVO PRO L EYE EVO PRO SANE PRO OUTDOOR	L a5'		PPTION 1 C	DPTION 2 	OPTION 1	OPTION 2		INFRAME GLAZING
a192 TYC/ 	L EYE EVO PRO L EYE EVO PRO SANE PRO OUTDOOR	S a5 <sup>5</sup>		a714		a731	a779	a540	
a193 EVIL	L EYE EVO PRO L EYE EVO PRO SANE PRO OUTDOOR	: L			a715				
	L EYE EVO PRO ::					a731	a779	a541	
	CANE PRO OUTDOOR	S		a715	a723	a731		a542	
a194 EVIL			a516	a715	a723	a731		a543	
a196 TYC <i>i</i>	PANE DDO OUTDOOD	L a5′		a714	a715	a731	a779	a540	
a197 TYC <i>i</i>	CANE PRO OUTDOOR	S a5′	15 / a516	a714	a715	a731	a779	a541	
a198 EVIL	L EYE HALFRIM PRO	S a5′	 15 / a516	a714	a715	a731	a779	a531	
a199 EVIL	L EYE HALFRIM PRO	S a5	 15 / a516	a708	a715	a731	a779	a539	
a402 EVIL	L EYE HALFRIM	L a51		a714	a715	a731	a779	a530	
a403 EVIL	L EYE HALFRIM	S a51	 15 / a516	a714	a715	a731	a779	a531	*
a412 EVIL	L EYE HALFRIM X	S a5	 15 / a516	a708	a715	a731	a779	a539	*
a418 EVIL	L EYE EVO		a516	a715	a723	a731		a542	*
a419 EVIL	L EYE EVO	S	a516	a715	a723	a731		a543	
a421 KUM	MACROSS HALFRIM		a515	a714	a708	a779		a544	
a423 WHI	IPSTART		a515	a714	a708	a779			•
a424 KUM	MACROSS 2.0		a515	a714	a708	a779			•
a425 WILI	DCHARGE		a515	a714	a708	a779			•
a426 MAT	TIC		a515	a714	a708	a779			•
a427 3MA	ATIC		a515	a714	a708	a779			•
a428 EXC	CALATE		a515	a714	a708	a779			•
a429 SPR	RUNG		a515	a714	a708	a779			•